

## Stylish, customizable, all-in-one ENT procedure system with Futura Evo Power Procedure Chair



The Euroclinic Otocompact Steel joins great functionality with a sophisticated look to significantly improve ENT workflow. Its modular construction and mobility replace built-in cabinetry and inefficiently designed work spaces. The Steel combines commonly used elements of an outpatient setting with cutting-edge technology for excellent performance.


The Steel integrates wireless controls, instrument tablet, UV disinfection storage, LED light sources, and wireless endoscopy video system into one ergonomic, attractive, and easy-to-use unit. With a top of the line ENT chair and precision-made mobile cabinet housing constructed from the highest quality materials, the Steel is truly a complete office solution that is excellently engineered. Whether you are starting a new clinic or replacing several items, the Steel is the ideal option for customization.



**Instrument Panel**



**Futura Evo  
Procedure Chair**

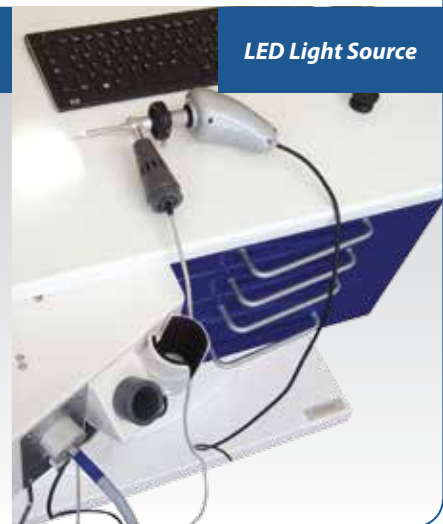

**LED Light Source**

### Key Components

#### Integrated Italian Leather Seating

- » Futura Evo ENT chair with four powerful motors, three electronic memory settings, wireless foot control, and non-deformable foam padding
- » SPD/A doctor's stool built from stainless steel with height adjustment and swiveling back rest

#### Mobile European Styled Workstation

- » Instrument management system with electronic console for dual control of chair and instruments
- » UV lamp germicide compartment kills common bacteria and viruses
- » Metal trays for instrument storage
- » Powerful suction system with dry pump
- » Insufflation system with nasal olive and autoclavable cannula
- » Auricular washing syringe with anti-sprinkle shield and microprocessor-controlled heating element
- » Larynx mirror heater programmable for 10-120 seconds
- » Air compression system for medications with 2 atomizers for liquids and powders
- » Antimicrobial filter is 99.99% effective against viruses and bacteria
- » Endoscope holders for your flexible and rigid endoscopes

#### Advanced Lighting System

- » Bright LED light source (180 Lumen at 5,500° K) with direct connection to rigid endoscopes
- » Neutral LED light source (180 Lumen at 5,000° K) with Storz connector

#### Integrated Video Endoscopy Suite

- » Wireless Camera EVS ED420 with ergonomic design and wireless foot pedal for hands-free photos and video
- » 23" HP TouchSmart PC monitor on articulating arm with wireless keyboard and mouse
- » Video capture and data management software
- » Reader/Writer smart card kit for easy record keeping



### Futura Evo Specifications

Power Supply Tension	110V or 230V-50/60Hz
Linak Motor Low Tension	24 Vdc
Motor Lifting Push	6000N fast run
Back Rest Motor Push	6000N
Seat Motor Push	6000N
Unique Leg Rest Motor Push	6000N
Lifting Capacity	180 kg
Absorbed Power	120 VA
Weight	120 kg
Base Length	800 mm
Maximum Height of Leg Rest	1250 mm



### Otocompact Steel Specifications

Power Supply Tension	230 V +/- 10% optional 110V +/-10%
Frequency	50/60Hz
Absorbed Power (MAX configuration)	1200 VA
Weight	100 kg
Max. Number of Instruments	4
Insufflator Air Capacity	26 l/min (variance until 250 Kpa)
Washing Syringe Capacity	350/400 cc/min (a 2 Bar)
Nominal Suction (motor)	-86 KPa (-650 mmHg) 50 L/min
Suction Jar	1 Lt (optional 2 Lt) multiple or single use
Power Supply Fuses	F1, F2=6, 3 AT